# *Avicenum* **ORTHO 360**

## **(EN) INSTRUCTION MANUAL**

Avicenum ORTHO 360 knee bandage Type 02

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Avicenum ORTHO 360 Kniebandage Typ 02

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Avicenum ORTHO 360 bandáž kolenná typ 02

**(EN) Dear customer,**

In your hands, you are now holding a knee bandage for post-operative and post-traumatic therapy. You can maximize the benefits of the product by its proper use and maintenance.

**Avicenum ORTHO 360 knee bandage Type 02** is manufactured by the circular knitting technology. The materials used are non-irritating and latex-free. The product achieves the desired therapeutic effect by exerting compression pressure in the knee area (you must choose the correct size). The therapeutic effect of the bandage is enhanced by the combination with a hook-and-loop fastener. It is intended for knee support in post-traumatic and post-operative therapy. The product is most effective during physical strain. The product is also suitable for sports activities. In the upper part, the product has a hem with silicone studs that prevent sliding. You should remove the product during longer breaks. **Avicenum ORTHO 360 knee bandage Type 02** is a medical device intended for use both at home and in medical facilities.

**Product name:**

**Avicenum ORTHO 360 knee bandage Type 02**

**Product variants:**

- fastening strap, size S–XXL

**User profile:**

**Avicenum ORTHO 360 knee bandage Type 02** is a medical device intended for use both at home and in medical facilities.

**Intended use:**

A medical device designed to support the knee during post-traumatic and post-operative therapy, in case of flabby muscles, inflammation, or injury to knee ligaments and tendons and to prevent swelling. This product is also suitable for sports activities as a therapeutic or preventive aid.

**Indications:**

- knee support in post-operative or post-traumatic therapy
- instability and weakness of knee ligaments
- sprain or injury of knee ligaments and tendons
- inflammation and swelling in the knee area
- prevention of injury and knee strengthening (e.g., during sporting activities)

If this is the first time you use a knee bandage, consult your doctor or have your medical condition regularly examined.

**Contraindications:**

- damaged skin
- acute weeping skin manifestations
- lymphatic swelling
- hypersensitivity to materials used

- limb sensitivity disorders (e.g., peripheral neuropathy)

Your doctor should consider the use of the product in view of the known contraindications and possible risks compared to the product's therapeutic benefits.

**Target group:**

Patients with knee instability, patients undergoing post-operative or post-traumatic therapy of the knee, and patients with the above indications. The product is intended for use on healthy, intact skin. The product is intended for all age groups and both sexes. The product may be used in pregnancy without limitation.

**Side effects:**

Side effects may occur especially when using a wrong size of the product or if the product is used outside strenuous activities. If the product is too tight or too tight support tapes, there may be pressure sores, impaired blood circulation due to compression of blood vessels, or deterioration of sensitivity due to excessive compression of the nerves. Sensitive and hypersensitive individuals may experience skin irritation (redness, itching, or blistering) at the point of contact of the product with the skin.

**How to choose the correct size:**


In order to achieve the desired effects, you must select the correct size of the product. Determine the size by measuring the circumference at the points indicated in the figure.

**Note:**

Any swelling or instability of the knee should be examined by your doctor. If this is the first time you experience swelling or instability of the knee, consult your doctor before using the product. Do not use the product if you do not know the cause of the swelling or instability of the knee. Always consult your doctor before combining the product with other medical devices (e.g., compression stockings).

- use a tape measure
- take the measurements standing up, preferably in the morning right after waking up
- get assistance from another person
- measured values may change during the therapy – always verify the correct size before using a new product

Size	Circumference below the knee in cm (point d)	Circumference above the knee in cm (point f)
S	30 - 34	41 - 46
M	34 - 38	46 - 51
L	38 - 42	51 - 56
XL	42 - 46	56 - 61
XXL	46 - 50	61 - 66

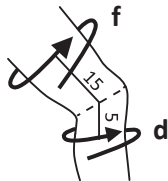**If you are on the border between two sizes**, choose the bigger one. If you have difficulties choosing the size of the product, visit a specialized medical supplies store or directly the manufacturer.

**If you have difficulties putting on the product**, visit a specialized medical supplies store and request a demonstration of the correct method of applying it.

**With a correctly selected size:**

- you should feel noticeable but not unpleasant pressure at the knee area
- the product does not slide or roll down
- the product does not choke at any point

Size and manufacturer information can be found on the tag sewn at the top edge of the product.

**Method of use:**

Put on the product by gradually pulling it over the knee. Never put on the bandage by merely pulling the top edge – always put the bandage on gradually. When putting on, the product may not come into contact with sharp objects (e.g., jewellery) or long nails. Stretch the product evenly without creases for effective distribution of pressure. To increase the therapeutic effect, put a support tape over the bandage. You can see the method of application in the picture:

If you use body cosmetics, wait about 15 minutes after application before you put on the product.

In the case of visible mechanical damage to the product, stop wearing it and preferably replace it with a new product.

**Care and maintenance:**

**Avicenum ORTHO 360 knee bandage Type 02** is a medical device intended for repeated use by a single person. The product will maintain its properties and safety for 12 months from the first use under these conditions:

- the product can be washed in a washing machine (maximum washing temperature 30 °C, gentle cycle) in a soapy solution or using a special preparation for washing compressing stockings and sleeves using a wash bag
- do not use fabric softeners
- do not use centrifugation; squeeze the wet product between two towels
- dry the product in a horizontal position away from a direct heat source (radiator, sun)

- do not iron the product
- the maximum number of washing cycles is 120
- store in a dry and dark place, preferably in the original packaging
- do not interfere mechanically in the product, do not repair places damaged by use, do not cut the edges of the product
- pay special attention to the hook-and-loop fastener on the support tape: contact of the hook-and-loop fastener with the knit area may cause damage
- contact of the hook-and-loop fastener with other materials (fabric, textiles, dirt and grime) may damage it and cause it to stop working in whole or in part
- the product may not come into contact with organic solvents or bleaching agents

The product is safe to use for 60 months from the date of manufacture. This date is included in the product's batch number in the format YY/MM-1234567 where YY is the year and MM the month of manufacture.

#### Care symbols:

Do not use fabric softeners!

#### Prohibited use:

Stop using the product immediately if:

- If it shows signs of mechanical damage (replace the product)
- Your health has deteriorated which could be associated with the use of the product – consult your doctor at once

#### Composition:

60 % PAD Nylon, 40 % elastan LYCRA®, silikon

This medical device does not contain an active substance, including human blood or plasma derivatives, tissues or cells of human origin or derivatives thereof, tissues or cells of animal origin or derivatives thereof.

#### Disposal:

The product can be disposed of with normal household waste or in a textile waste container. Although the product is intended for use on healthy and intact skin, it may be contaminated with body fluids during use. In that case, dispose of the product as an infectious material.

The material used is not toxic unless burned. Do not burn the used bandage.

**Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.**

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**Manufacturer:**

ARIES, a.s., Studenec 309, 512 33 Studenec, Czech Republic  
Last reviewed: 17.07.2020

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